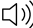
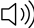
# Questions & Answers

## 1- Does "LO" means "Low Battery"?

**No**, it means the temperature is below 34°C (93.2°F) in forehead mode. Simply bring the thermometer closer to the forehead: 3-5cm (1.2-2in)

## 2- How to turn the sound On/Off?

1. When the device is On, press the "mem" key to set the sound to On or Off.
2. The LCD screen will show  and you will hear a short beep.
3. To turn the sound Off, press "mem" button again,  will change into .

## 3- How to switch from °C to °F

When the device is Off, press the "mem" button for six seconds. You can now switch between Fahrenheit degree (°F) and Celsius degree (°C).

## 4- How to go to the Object mode?

1. Press and hold  for 6 seconds.
2. The forehead mode icon  will disappear, and the object mode icon will appear .

## 5- Why does the screen displays "HI"?

It means that the target temperature is higher than the maximum of the measurement range: 43°C in Forehead Mode and 93.2°C in Object Mode.

## 6- The results are different compared to other thermometers. Is it normal?

**YES.** It is normal to obtain different results depending on the type of thermometer used (frontal, ear, oral or rectal).

**PLEASE REFRAIN FROM COMPARING RESULTS** and instead refer to the user's manual.

## More Questions?

Contact us at [Support@ByFloProducts.com](mailto:Support@ByFloProducts.com)

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