



Quick Start

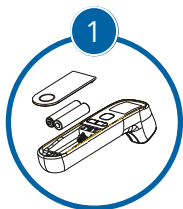



Unpack thermometer and insert batteries.

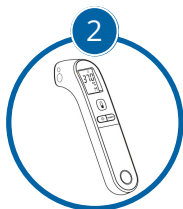

Place the thermometer at room temperature for 30 minutes before using.

FOREHEAD MODE

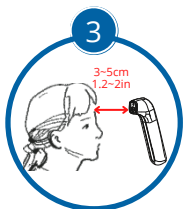

Point the thermometer 3-5cm away from the forehead and press .



Normal

34~37.1°C
93.2~98.8°F



Attention

37.2~38.1°C
98.9~100.6°F




Fever

38.2~43°C
100.7~109.4°F

OBJECT MODE



Press and hold  for 6 sec.

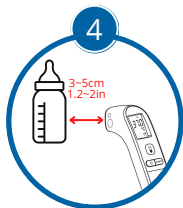

Point the thermometer 3-5cm away from the object and press .




CAUTION!


It is **normal to obtain different results** depending on the type of thermometer used (frontal, ear, oral or rectal).

Therefore, **THE READINGS ON THE FOREHEAD WILL NOT BE THE SAME AS ORAL OR AXILLARY READINGS.**

Démarrage Rapide

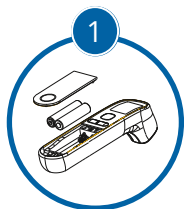


Déballiez le thermomètre et insérez les piles.

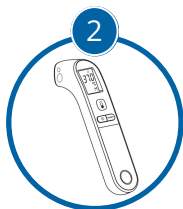

Placez le thermomètre à température ambiante pendant 30 minutes avant utilisation.

MODE FRONTAL

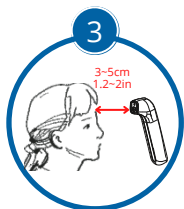

Pointez le thermomètre à 3-5cm du front et appuyez sur .



Normal

34~37.1°C
93.2~98.8°F



Attention


37.2~38.1°C
98.9~100.6°F




Fever


38.2~43°C
100.7~109.4°F

MODE OBJET

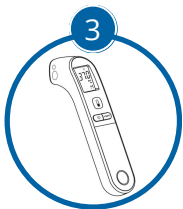

Appuyez et maintenez  pendant 6 sec.

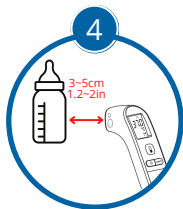

Pointez le thermomètre à 3-5cm de l'objet et appuyez sur .



ATTENTION!

Il est normal d'obtenir des résultats différents selon le type de thermomètre utilisé (frontal, auriculaire, oral ou rectal).

Par conséquent, **LES LECTURES FRONTALES NE SERONT PAS LES MEMES QUE LES LECTURES ORALES OU AXILLAIRES.**